# ***virtufit***

**Wireless Neckmassager Pro**

VF11009

**USER MANUAL**

📷 @virtufit\_fitness

📺 VirtuFit

🌐 [www.virtufit.com](http://www.virtufit.com)


## HOW TO USE

**NECK**


**BACK**

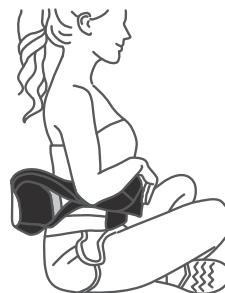
**WAIST**

**ABDOMEN**

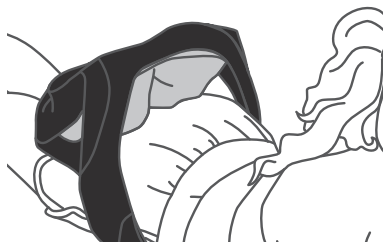**HIPS**

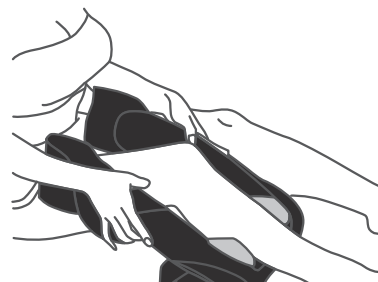**CALF**

**FOOT**

**ARMS**

**HIGH INTENSITY**

## SAFETY INSTRUCTIONS

### ⚠ WARNING!

*In case of problems with your body or skin, consult your doctor and follow his instructions. Do not use this device while sleeping. Do not use on damaged skin or damaged sensitive nerves. Always start with the lowest level and gradually increase the level.*

**NOTE!** Due to safety concerns, the product should not be used while driving.

Persons with the following condition should consult a doctor before using this device.

- Acute illness
- Malignant tumour
- Infectious disease
- Pregnancy
- Cardiac dysfunction
- High fever
- Abnormal blood pressure
- Skin sensitivity disorders or skin problems
- Receipt of medical treatment, especially for discomfortak

**NOTE!** Should not be used by persons in the first trimester of pregnancy, with a pacemaker or AICD, or who are being treated for an existing deep vein thrombosis.

change the direction of rotation of the massager to change and maintain.

### Speed

3 different speeds from low to high, choose the suitable speed according to your need.


### Heat on/off

Press this button to turn the heat function on or off while the massage continues.

## CONDITION OF USE

- Store in a dry and dust-free place
- Do not spill any liquid on the device
- Do not immerse the appliance in water
- Do not clean with chemicals

## USE
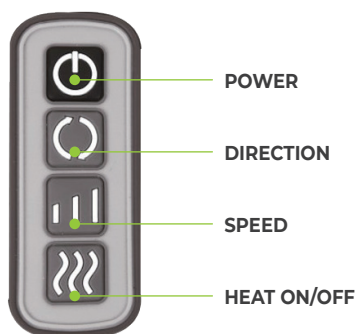
### Power

- Press to switch the massager on.
- The automatic timer switches the device off after 15 minutes.
- Press the on/off button at any time to switch the appliance off.
- Heat will be automatically turned on.

### Direction

When the device is switched on, the massage roller automatically changes its rolling direction every 60 seconds. Press the button to

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